

School balance KERN EMS



Entry level model in the low-cost range with large weighing plate

Features

- Especially suitable for use in schools and universities, for example for biology, chemistry, physics
- Large, shock proof weighing plate made of plastic, with conductive lacquer
- Particularly flat design
- Ergonomically optimised key pad with large keys and a high-contrast LCD display
- Secure and non-slip positioning through rubber feet
- Adjusting program CAL for quick setting of the balance accuracy, external test weights at an additional price, see *Test weights*
- **1** Draught shield standard for models with weighing plate size **A**, weighing space W×D×H 145×145×65 mm
- Suitable for common school LIMS systems


Technical data

- LCD display, digit height 25 mm
- Dimensions of weighing surface, **A** Ø 105 mm **B** W×D 175×190 mm, see larger picture
- Weighing plate material **A** plastic, with conductive lacquer **B** plastic
- Overall dimensions W×D×H 200×280×63 mm (without draught shield)
- Optional battery operation, 9 V block not included in scope of delivery, operating time up to 40 h
- External mains adapter standard
- Net weight approx. 1,4 kg
- Permissible ambient temperature 5 °C/35 °C

Accessories

- **2** Stainless steel weighing plate, only for models with weighing plate size **B**, KERN EMS-A01

STANDARD

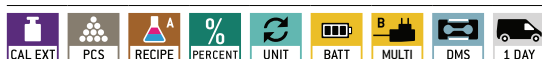

OPTION



Model	Weighing capacity [Max]	Readability [d]	Reproducibility	Linearity	Weighing plate	Option	
						DAKKS Calibr. Certificate	
KERN							
EMS 300-3	300	0,001	0,002	± 0,005	A	963-127	
EMS 3000-2	3000	0,01	0,02	± 0,05	B	963-127	
EMS 6K0.1	6000	0,1	0,1	± 0,3	B	963-128	
EMS 12K0.1	12000	0,1	0,1	± 0,3	B	963-128	
EMS 6K1	6000	1	1	± 3	B	963-128	
EMS 12K1	12000	1	1	± 3	B	963-128	

Pictograms

	Internal adjusting: Quick setting up of the balance's accuracy with internal adjusting weight (motordriven)		KERN Communication Protocol (KCP): It is a standardized interface command set for KERN balances and other instruments, which allows retrieving and controlling all relevant parameters and functions of the device. KERN devices featuring KCP are thus easily integrated with computers, industrial controllers and other digital systems		Suspended weighing: Load support with hook on the underside of the balance
	Adjusting program CAL: For quick setting up of the balance's accuracy. External adjusting weight required				Battery operation: Ready for battery operation. The battery type is specified for each device
	Easy Touch: Suitable for the connection, data transmission and control through PC or tablet.		GLP/ISO log: The balance displays serial number, user ID, weight, date and time, regardless of a printer connection		Rechargeable battery pack: Rechargeable set
	Memory: Balance memory capacity, e.g. for article data, weighing data, tare weights, PLU etc.		GLP/ISO log: With weight, date and time. Only with KERN printers.		Universal plug-in power supply: with universal input and optional input socket adapters for A) EU, CH, GB; B) EU, CH, GB, USA; C) EU, CH, GB, USA, AUS
	Alibi memory: Secure, electronic archiving of weighing results, complying with the 2014/31/EU standard.		Piece counting: Reference quantities selectable. Display can be switched from piece to weight		Plug-in power supply: 230V/50Hz in standard version for EU, CH. On request GB, USA or AUS version available
	Data interface RS-232: To connect the balance to a printer, PC or network		Recipe level A: The weights of the recipe ingredients can be added together and the total weight of the recipe can be printed out		Integrated power supply unit: Integrated in balance. 230V/50Hz standard EU. More standards e.g. GB, USA or AUS on request
	RS-485 data interface: To connect the balance to a printer, PC or other peripherals. Suitable for datatransfer over large distances. Network in bus topology is possible		Recipe level B: Internal memory for complete recipes with name and target value of the recipe ingredients. User guidance through display		Weighing principle: Strain gauges: Electrical resistor on an elastic deforming body
	USB data interface: To connect the balance to a printer, PC or other peripherals				Weighing principle: Tuning fork: A resonating body is electromagnetically excited, causing it to oscillate
	Bluetooth* data interface: To transfer data from the balance to a printer, PC or other peripherals		Totalising level A: The weights of similar items can be added together and the total can be printed out		Weighing principle: Electromagnetic force compensation: Coil inside a permanent magnet. For the most accurate weighings
	WiFi data interface: To transfer data from the balance to a printer, PC or other peripherals		Percentage determination: Determining the deviation in % from the target value (100 %)		Weighing principle: Single cell technology: Advanced version of the force compensation principle with the highest level of precision
	Control outputs (optocoupler, digital I/O): To connect relays, signal lamps, valves, etc.		Weighing units: Can be switched to e.g. nonmetric units. See balance model. Please refer to KERN's website for more details		Verification possible: The time required for verification is specified in the pictogram
	Analogue interface: to connect a suitable peripheral device for analogue processing of the measurements		Weighing with tolerance range: (Checkweighing) Upper and lower limiting can be programmed individually, e.g. for sorting and dosing. The process is supported by an audible or visual signal, see the relevant model		DAKkS calibration possible (DKD): The time required for DAKkS calibration is shown in days in the pictogram
	Interface for second balance: For direct connection of a second balance		Hold function: (Animal weighing program) When the weighing conditions are unstable, a stable weight is calculated as an average value		Factory calibration (ISO): The time required for Factory calibration is shown in days in the pictogram
	Network interface: For connecting the scale to an Ethernet network		Protection against dust and water splashes IPxx: The type of protection is shown in the pictogram.		Package shipment: The time required for internal shipping preparations is shown in days in the pictogram
					Pallet shipment: The time required for internal shipping preparations is shown in days in the pictogram

*The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by KERN & SOHN GmbH is under license. Other trademarks and trade names are those of their respective owners.

KERN – Precision is our business

To ensure the high precision of your balance KERN offers you the the appropriate test weight in the international OIML error limit classes E1-M3 from 1 mg – 2500 kg. In combination with a DAKkS calibration certificate the best pre-requisite for proper balance calibration.

The KERN DAKkS calibration laboratory today is one of the most modern and best-equipped DAKkS calibration laboratories for balances, test weights and force-measurement in Europe.

Thanks to the high level of automation, we can carry out DAKkS calibration of balances, test weights and force-measuring devices 24 hours a day, 7 days a week.

Range of services:

- DAKkS calibration of balances with a maximum load of up to 50 t
- DAKkS calibration of weights in the range of 1 mg – 2500 kg
- Volume determination and measuring of magnetic susceptibility (magnetic characteristics) for test weights
- Database supported management of checking equipment and reminder service
- Calibration of force-measuring devices
- DAKkS calibration certificates in the following languages DE, EN, FR, IT, ES, NL, PL
- Conformity evaluation and reverification of balances and test weights

Your KERN specialist dealer: